


**MAC-2014     Mobile Active Count**

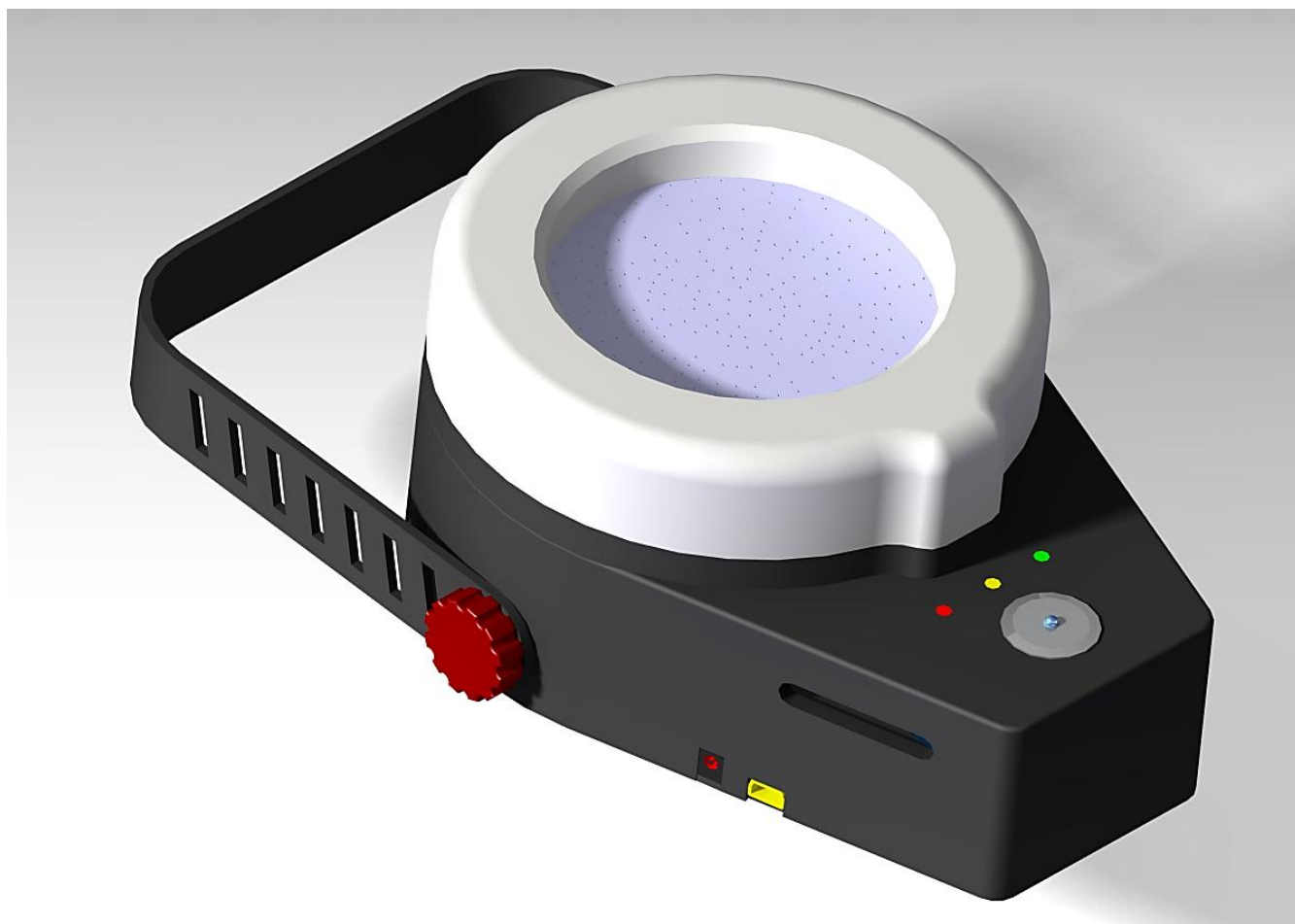

# **User Manual**

---

SISTEMA-MK GMBH

# User Manual MAC-2014

---

Sistema-MK GmbH

Ricarda-Huch-Str. 32

D-73760 Ostfildern

Phone / Fax: +49 711 57641776 / 50471652

E-Mail: [info \(at\) sistema-mk.eu](mailto:info@systema-mk.eu)

---

# Table of Contents

<b>MAC-2014 MOBILE ACTIVE COUNT</b>	<b>0</b>
<b>SYSTEM CONFIGURATION</b>	<b>4</b>
<b>MAC-2014 OPTIONS</b>	<b>4</b>
MAC-2014 STANDARD	5
Dimension:	5
MAC-2014 EXTENDED VERSION	6
Dimension:	6
<b>QUICK START CONFIGURATION</b>	<b>7</b>
<b>TRAFFIC LIGHT</b>	<b>7</b>
<b>FIRST STEP TO START UP</b>	<b>8</b>
Start MAC 2014	8
Start Printer	8
Start tablet	8
<b>SAMPLING MODE DEFINITION</b>	<b>9</b>
Auto start Activation	9
Print Text Activation	10
Print Text Configuration	11
Delay time	11
Location	11
Measurement volume	12
Measurement Start	12
<b>ENVIRONMENT MEASUREMENT</b>	<b>13</b>
<b>CO<sup>2</sup> MEASUREMENT</b>	<b>14</b>
<b>MEASUREMENT DATA LOG FILE</b>	<b>15</b>
<b>PRINTER SETUP</b>	<b>16</b>
<b>CALIBRATION SETUP</b>	<b>18</b>

---

<b>NETWORK SETUP</b>	<b>20</b>
----------------------	-----------

---

HOT SPOT Setup	20
----------------	----

WLAN Setup	21
------------	----

## System Configuration


*MAC-2014 is available in two versions, Standard and incl. range extender.*

### MAC-2014 Options

	MAC-2014 Standard	MAC-2014 Range Extender
Battery capacity	1000mAh	6.800mAh
Auto Start Function	Standard	Standard
Altitude Measurements	Standard	Standard
Atmospheric Pressure	Standard	Standard
Humidity Measurements	Optional	Standard
Temperature Measurements	Optional	Standard
CO <sup>2</sup> Measurement	Optional	Optional
Printer Functionality	Optional	Optional
Real Time Clock	Optional	Optional
Camera Stand Connection Threat	Optional	Standard
HOT Spot	Standard	Standard
WLAN Function	Optional	Optional



## SYSTEM CONFIGURATION

### MAC-2014 Standard

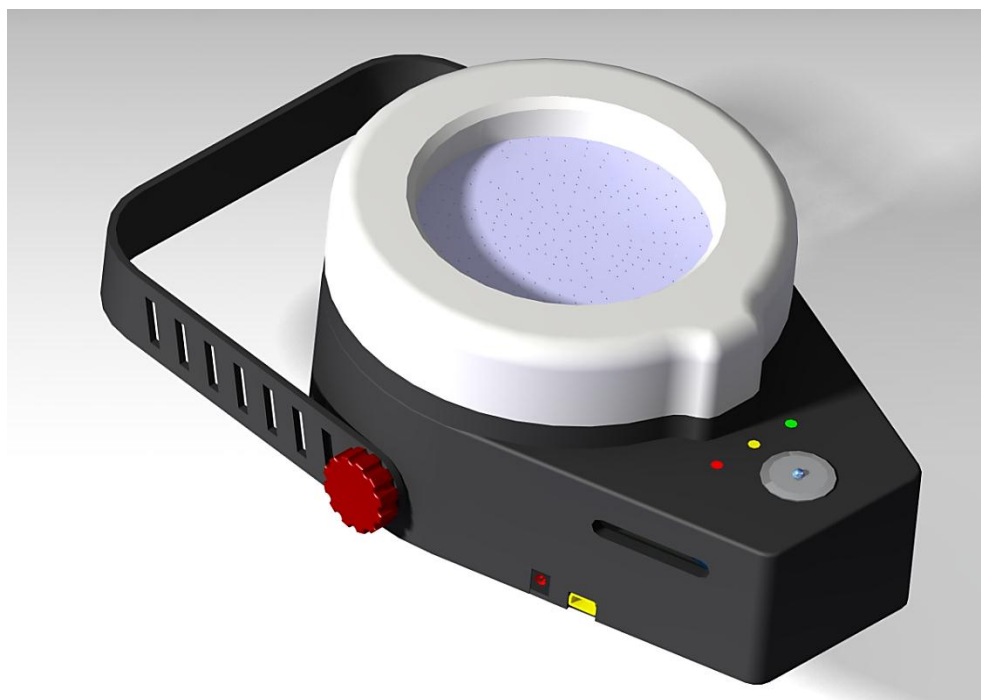

MAC-2014 as Standard version


DIMENSION:

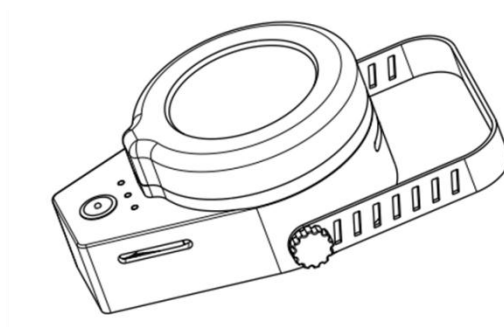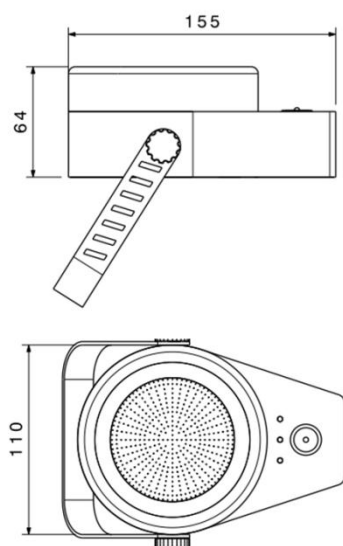

MAC 2014 Standard



## SYSTEM CONFIGURATION

### MAC-2014 Extended Version

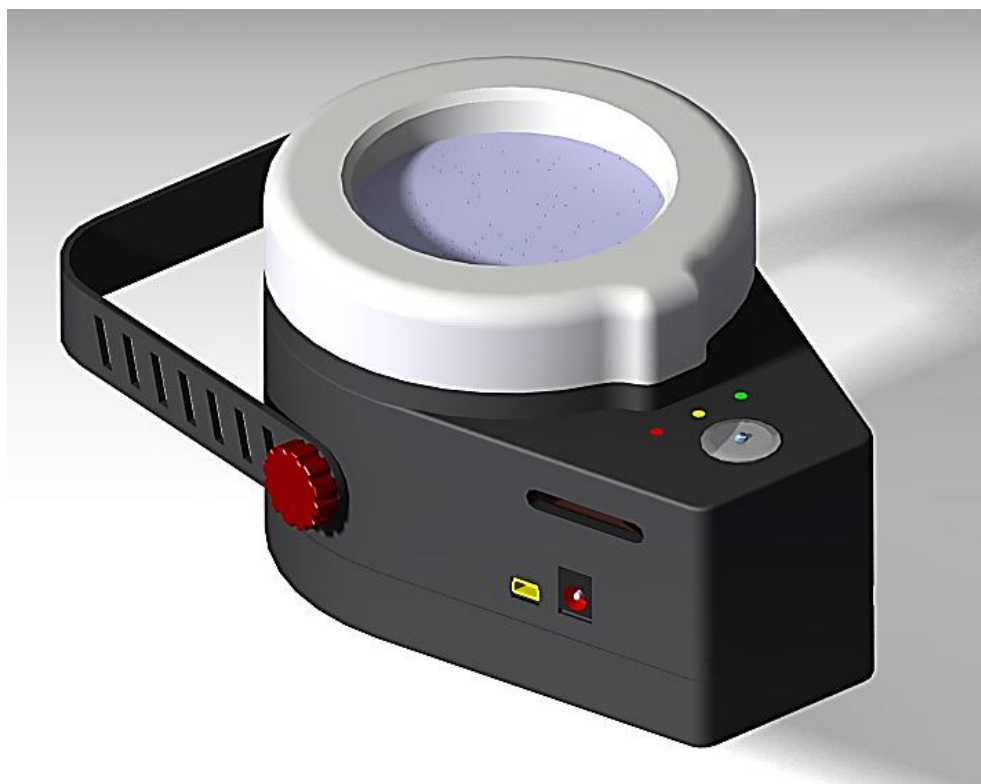

MAC-2014 as Expended version

DIMENSION:

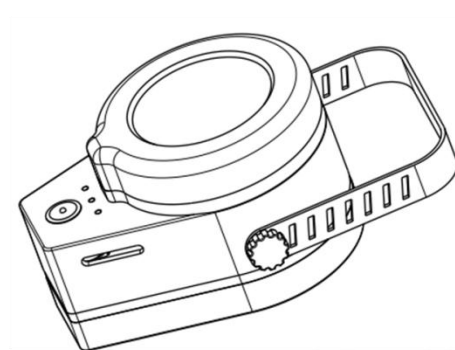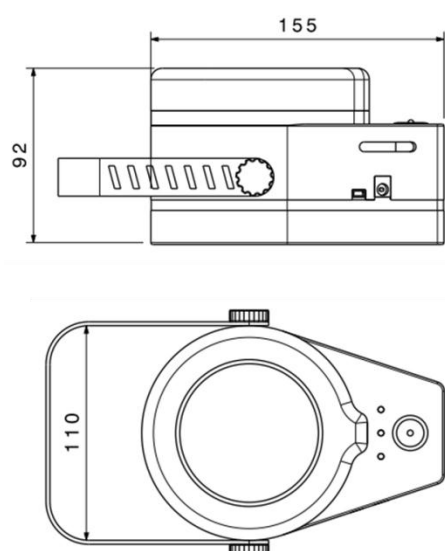
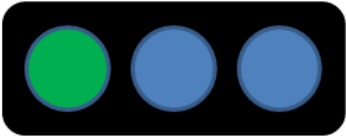
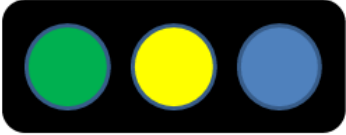
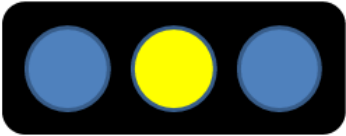




MAC 2014 Range Extender

## Quick Start Configuration

### Traffic Light

MAC-2014 is working with 3 status light green, yellow and red.

Status Light	Definition
	Self-Test Status All lights on –everything is OK
	System read System is ready to measure
	System waiting time System has started – waiting time
	Sampling mode System is currently sampling
	Sampling mode finished Sampling is finished please remove plate – yellow blinking
	System failure – red blinking



## QUICK START CONFIGURATION

### First Step to Start Up

#### START MAC 2014

Push the blue button of the MAC 2014, as soon as the green light is on the MAC is ready to work.

#### START PRINTER

Push the start button of the printer the printer will automatically connect to the MAC 2014 – if the blue wireless sign is on the MAC 2014 is successfully connected.

#### START TABLET

Push the start button of the tablet PC.

Connect the WLAN to the MAC\_2014 Hot Spot



Start the FireFox WEB Browser



Tip in the following network IP (192.168.42.1) in case the MAC was setup as Hotspot.

.

---

## Sampling Mode Definition

MAC\_2014, s/n: 00001  
[System Info](#) | [Print](#) | [Shutdown](#) | [Battery](#)

**LIGHTHOUSE**  
WORLDWIDE SOLUTIONS BENELUX BV

Start Measurement  
C° and H% Measurement  
CO<sup>2</sup> Measurement  
DataLog  
Printer Setup  
Calibration Setup  
Network Setup

☐ **Autostart Activation** (Measurement will automatic start after restart)

☐ **Print Text Activation**

Printtext:

ID: 151	Date: 24-08-2014	Time: 15:16:28
Volume: 50 Liter	Location:	

Delay Time :  Seconds  
Measurement :  Liters

Location: 1

Location: 2

Location: 3

**Measurement Start - Volume Definition -**

50	100	250	500	1000
----	-----	-----	-----	------

Measurement START

### AUTO START ACTIVATION


☐ **Autostart Activation** (Measurement will automatic start after restart)

System will start automatically after re start with the setting of the last measurement.

- ☐ Auto start Activation – OFF
- ☒ Auto start Activation – ON


## SYSTEM CONFIGURATION

### PRINT TEXT ACTIVATION

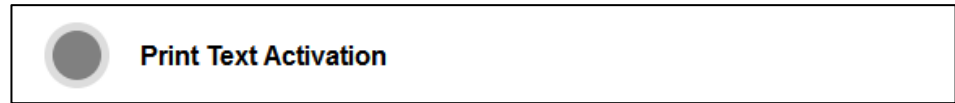

System will create a label text print file. The label will be printed as soon as the printer is switched on and connected by WLAN..

- ☐ Print Text Activation – OFF
- ☒ Print Text Activation – ON

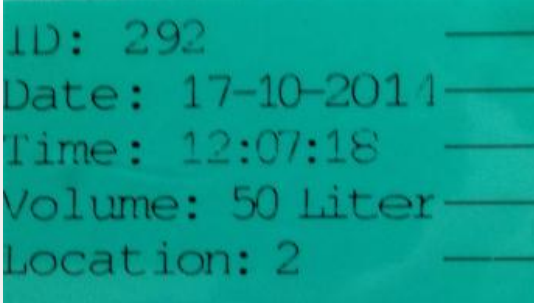A rectangular label with a light blue background and black text. The text is arranged in five lines, each followed by a horizontal line. The text reads: 'ID: 292', 'Date: 17-10-2014', 'Time: 12:07:18', 'Volume: 50 Liter', and 'Location: 2'.

```
ID: 292 _____
Date: 17-10-2014 _____
Time: 12:07:18 _____
Volume: 50 Liter _____
Location: 2 _____
```

## QUICK START CONFIGURATION

### PRINT TEXT CONFIGURATION

<b>Printtext:</b>	ID: 151    Date: 24-08-2014    Time: 15:16:28		
	Volume: 50 Liter	Location:	

ID:                      Counting number of Measurements

Date:                    Date of sampling

Time:                    Time of sampling

Volume:                Sampling Volume

Location:              Location of sampling

Print text will be automatically written into the print text field, in case of any needed of modification, the text can be changed by single activation of the text field.

### DELAY TIME

<b>Delay Time :</b>	<input type="text" value="1"/>
---------------------	--------------------------------

Definition of the delay time before sampling mode will start

### LOCATION

Location: 1	▲ ▼
Location: 2	
Location: 3	

Definition of the measurement location, any customer specific location can be typed in manual into the print text file.

If a customer specific location is required please type in the location name into the location field, the new location name will be shown in the future as a kin of pre selection.

<b>Location:</b>
------------------

Location: Dusche
Location: Wohnzimmer
Location: 1
Location: 3
Location: 5
Location: 10
Location: 14
Location: 2
Location: 2
Location: 20
Location:

## QUICK START CONFIGURATION

### MEASUREMENT VOLUME

<b>Measurement :</b>	<b>50</b>
----------------------	-----------

The measurement volume can be defined by typing in the required volume into the Measurement volume field or by selection of pre-defined the measurement volumes.

Measurement Start - Volume Definition -				
50	100	250	500	1000

### MEASUREMENT START


<b>Measurement START</b>
--------------------------

Start of the measurement

During the measurement is running the screen is showing the red logo.



As soon as the measurement is finished the system is showing the green logo.

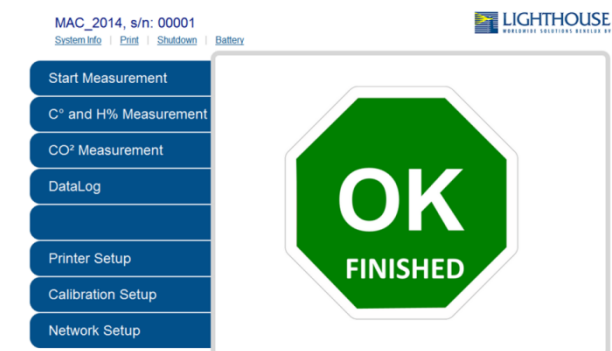

## Environment Measurement

MAC\_2014, s/n: 00001

[System Info](#) | [Print](#) | [Shutdown](#) | [Battery](#)

Start Measurement

C° and H% Measurement


CO<sup>2</sup> Measurement

DataLog

Printer Setup


Calibration Setup

Network Setup

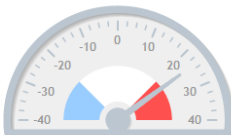

Temperature C

Air Temperature: 24 C

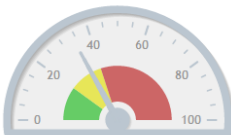

Air Humidity %


Air Humidity: 34 %

Barometric Pressure: 987 hPa      Altitude: 201 m

Air Temperature and Humidity Measurement

The measurement will take 10 seconds

Activation of the button “Air Temperature and Humidity Measurement” will start the environment measurement after 10 seconds blow time.

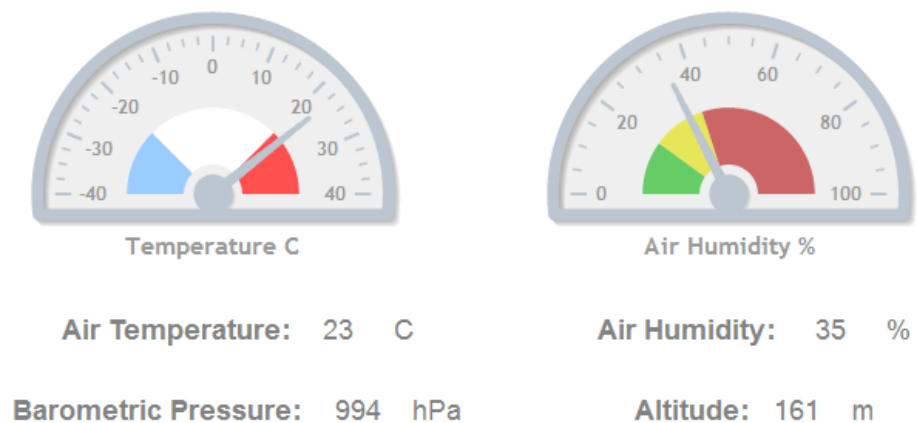

## CO<sup>2</sup> Measurement

MAC\_2014, s/n: 00001

[System Info](#) | [Print](#) | [Shutdown](#) | [Battery](#)

Start Measurement

C° and H% Measurement


CO<sup>2</sup> Measurement

DataLog

Printer Setup

Calibration Setup

Network Setup

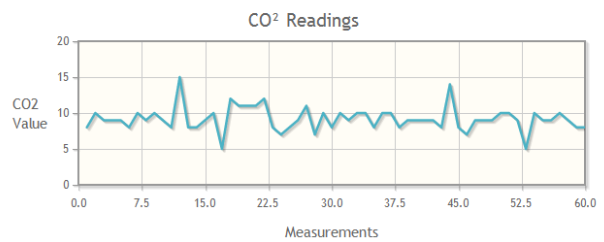CO<sup>2</sup> Measurement

The measurement will take 150 seconds  
90 seconds pre heating including blower  
60 seconds measurement without blower

## Measurement Data Log File

MAC\_2014, s/n: 00001

[System Info](#) | [Print](#) | [Shutdown](#) | [Battery](#)

Start Measurement

C° and H% Measurement

CO<sup>2</sup> Measurement

DataLog

Printer Setup

Calibration Setup

Network Setup


ID: 124	Date: 18-08-2014	Time: 14:34:34
Volume: 100 Liter	Location: 5	
ID: 125	Date: 18-08-2014	Time: 16:16:42
Volume: 50 Liter	Location: 2	
ID: 126	Date: 18-08-2014	Time: 16:22:47
Volume: 50 Liter	Location: 2	
Humidity: 34 % Temperature: 24 C Altitude: 196 m Barometric Pressure: 990 hPa		
ID: 127	Date: 19-08-2014	Time: 09:44:28
Volume: 100 Liter	Location: 2	
ID: 128	Date: 19-08-2014	Time: 10:34:55
Volume: 100 Liter	Location: 2	
ID: 128	Date: 19-08-2014	Time: 10:37:29
Volume: 50 Liter	Location: 2	
ID: 130	Date: 19-08-2014	Time: 10:43:44
Volume: 50 Liter	Location:	
ID: 130	Date: 19-08-2014	Time: 10:43:44
Volume: 50 Liter	Location:	
ID: 132	Date: 19-08-2014	Time: 10:44:17
Volume: 50 Liter	Location: 3	
ID: 132	Date: 19-08-2014	Time: 10:44:17
Volume: 50 Liter	Location: 3	
ID: 132	Date: 19-08-2014	Time: 10:44:17
Volume: 50 Liter	Location: 3	
ID: 135	Date: 19-08-2014	Time: 10:44:51
Volume: 50 Liter	Location:	
ID: 135	Date: 19-08-2014	Time: 10:44:51
Volume: 50 Liter	Location:	
ID: 137	Date: 19-08-2014	Time: 10:46:30



## Printer Setup

MAC\_2014, s/n: 00001

[System Info](#) | [Print](#) | [Shutdown](#) | [Battery](#)

Start Measurement	
C° and H% Measurement	
CO <sup>2</sup> Measurement	
DataLog	
Printer Setup	
Calibration Setup	
Network Setup	<div>Print Test String</div> <div>Printer Setup</div> <div>Network IP <input type="text" value="192.168.42.1"/></div>

Print label Sample:

```
ID: 292
Date: 17-10-2014
Time: 12:07:18
Volume: 50 Liter
Location: 2
```


Print successfully connected:



## SYSTEM CONFIGURATION

### PRINT LABEL TEXT

Free text label printing.

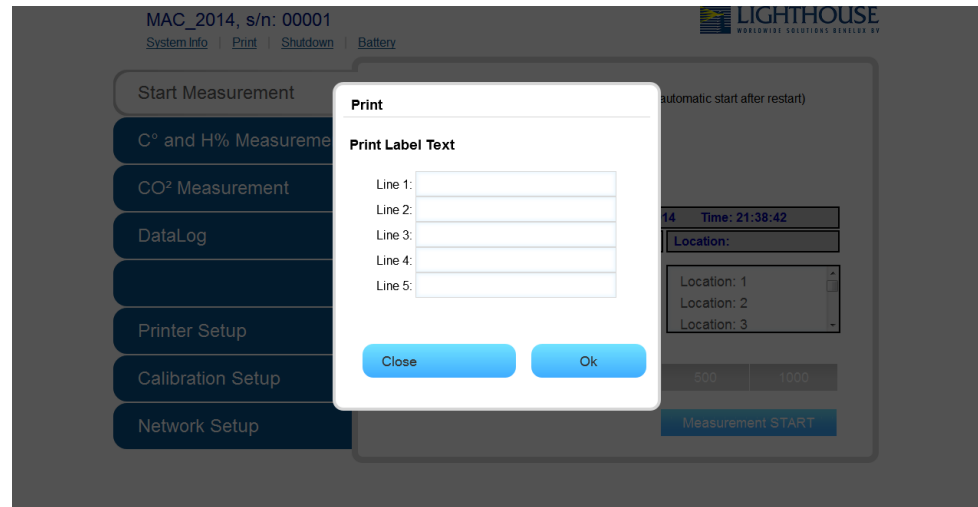

## Calibration Setup

MAC\_2014, s/n: 00001

[System Info](#) | [Print](#) | [Shutdown](#) | [Battery](#)

Start Measurement

C° and H% Measurement

CO<sup>2</sup> Measurement

DataLog

Printer Setup

Calibration Setup


Network Setup

P Value	300
I Value	1
D Value	2


Blower Level  Duration  Pressure Target

START Testrun

SAVE Setting

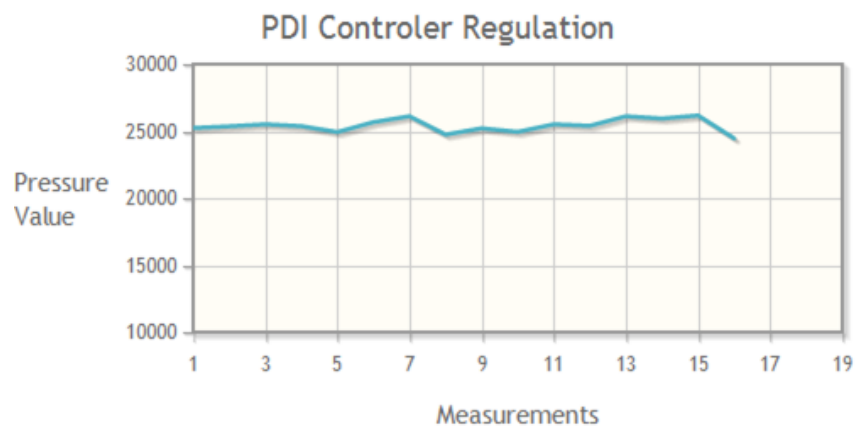

## SYSTEM CONFIGURATION

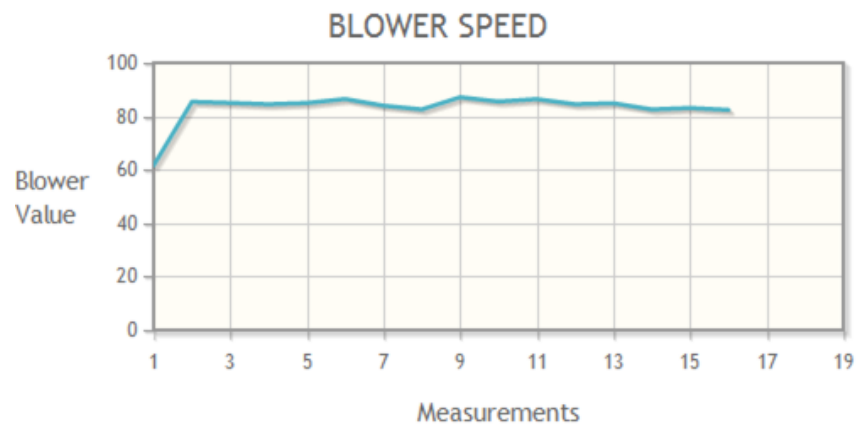



## Network Setup

MAC\_2014, s/n: 00001

[System Info](#) | [Print](#) | [Shutdown](#) | [Battery](#)



Start Measurement	<b>Hot Spot Setup</b>  Hotspot <b>WLAN Setup</b>  WLAN	<b>Aktiviere HotSpot</b>
C° and H% Measurement		<b>Aktiviere WLAN</b>
CO <sup>2</sup> Measurement		WLAN <input type="text" value="MAC_WEB"/> Password <input type="password" value="....."/>
DataLog		
Printer Setup		
Calibration Setup		
Network Setup		

### HOT SPOT SETUP



## SYSTEM CONFIGURATION

### WLAN SETUP



#### Remark:


In case the system was activated for WLAN Setup and the MAC is not close to any defined WLAN router, the system can always connect to the Tablet by activation the Mobile Hotspot function at the tablet. (Einstellungen, Mehr.. and Mobile Hotspot).

Select WLAN Button



## SYSTEM CONFIGURATION

Activate HotSpot

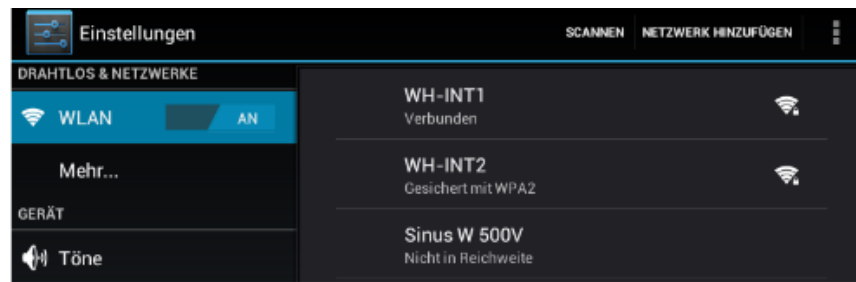

Login Information:

Netzwerk SSID: MAC\_2014

WPA2PSK: sistema-mk